



Long Implants


Product Catalog





Long Implants

Implant options up to 25mm in length

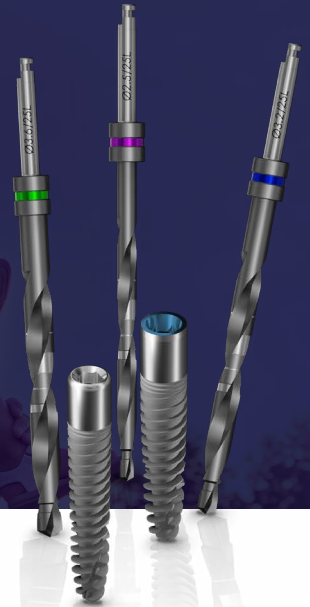


Adin's long implant range serves as an extension to the current Touareg-OS™ and CloseFit WP™ platform implants, allowing implant placement options up to 25mm in length.


In cases where more bone-to-implant contact is required, long implants may serve as a solution where the space allows for longer implants, alternatively, long implants can be used effectively as bi-cortical anchorage to compensate for lack of bone while still achieving implant stability in cases where implants of regular length might not.


Adin's long implant range carry specific features designed to enhance their performance such as narrow apex for limited spaces and a smooth neck to better suite gingival engagement.

Adin's unique and reliable OsseoFix™ surface treatment is utilized for both current Touareg™-OS and CloseFit™ WP platform long implants compatible drill set and depth probe have been designed to maximize ease of use and safety.

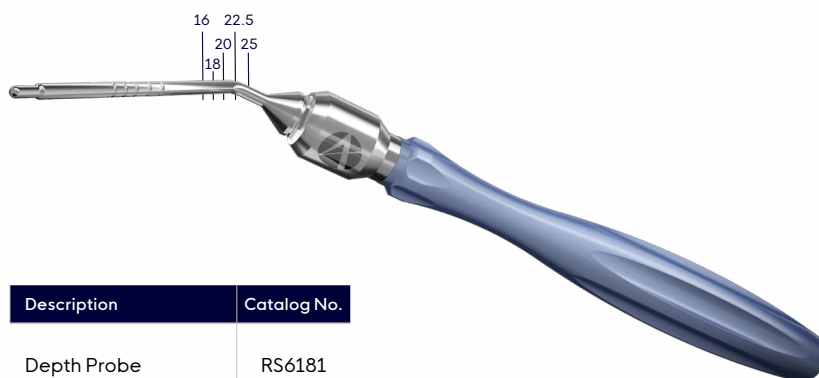


Description	Catalog No.
Depth Probe	RS6181
Surgical Driver	RS7001

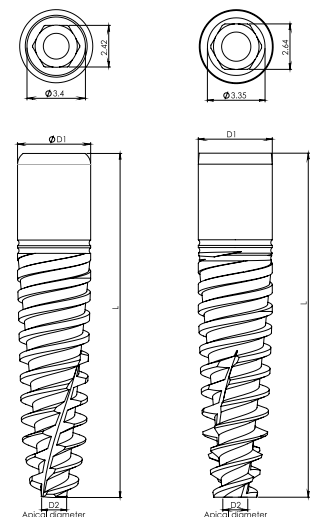

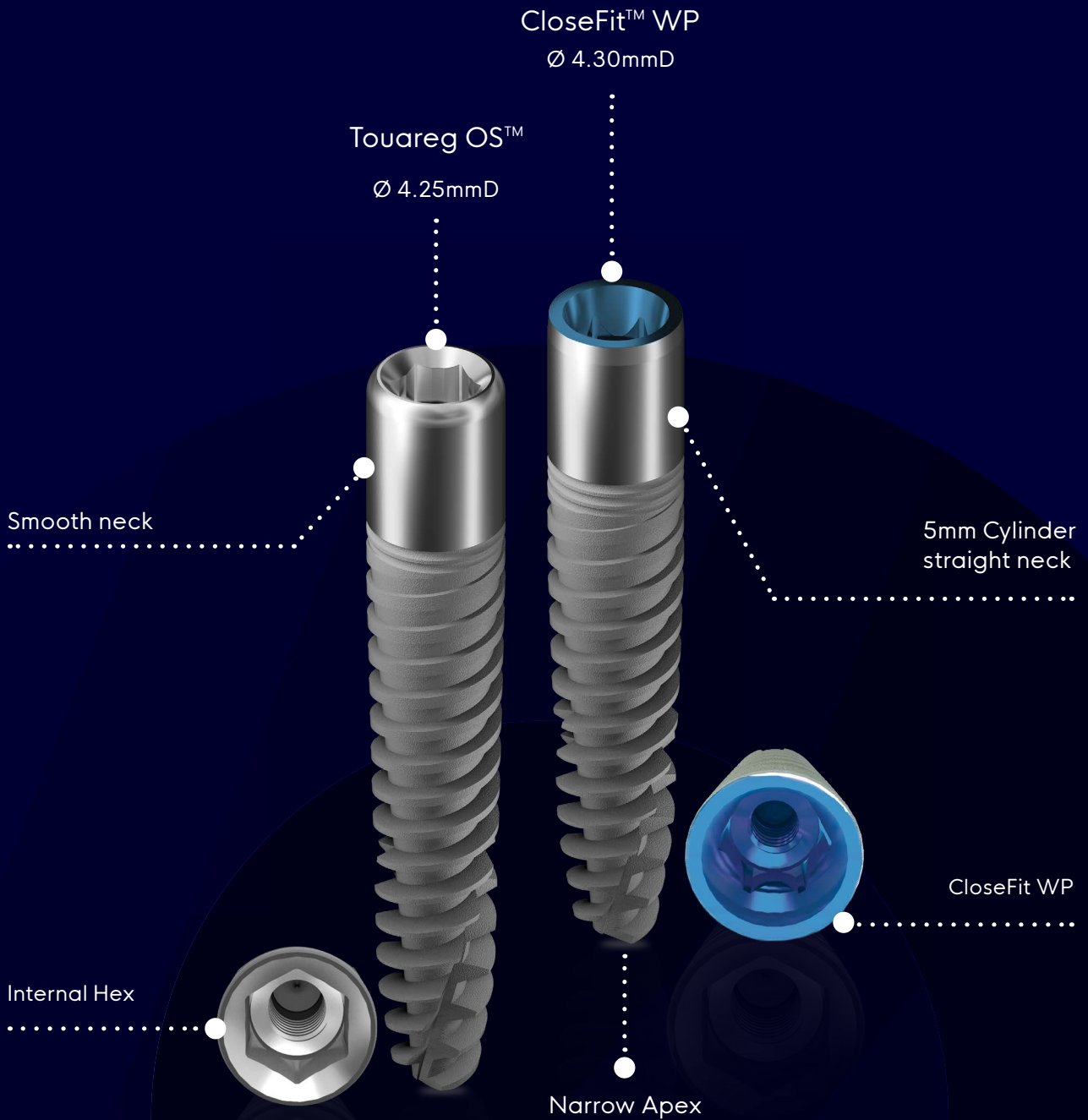


Touareg™- OS

Recommended insertion torque up to 50 Ncm


	D1	L	D2	Catalog No.
	4.2	20.0	1.5	ISF2042
		22.5	1.5	ISF2242
		25.0	1.5	ISF2542

CloseFit™ WP


Recommended insertion torque up to 50 Ncm

	D1	L	D2	Catalog No.
	4.3	20.0	1.5	ISPX2043
		22.5	1.5	ISPX2243
		25.0	1.5	ISPX2543

Touareg™-OS

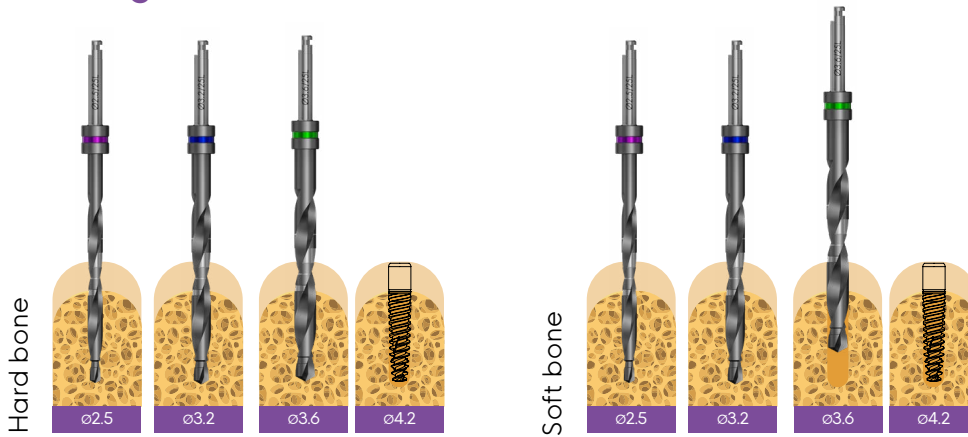

CloseFit™ WP

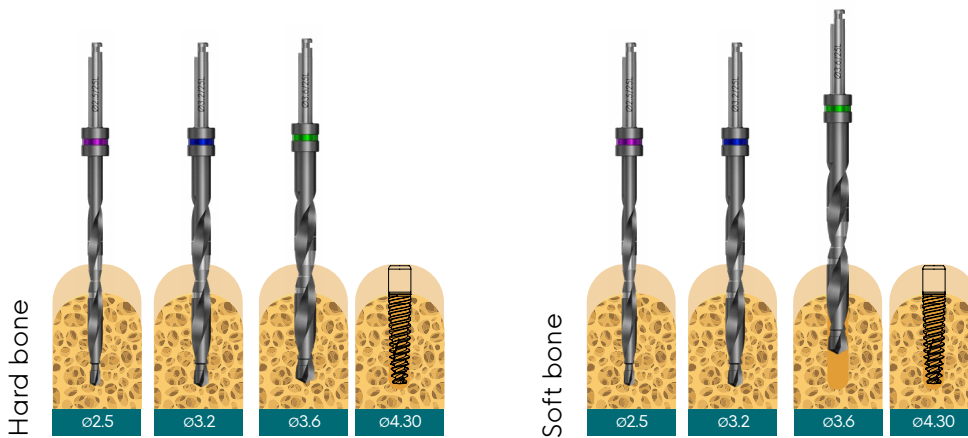

All measurements are in mm.

*For initial drilling, you may use Tri-step drill instead of ø2, ø2.8 and ø3.2 drills in sequence.

Note: (x.xx) Drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant. * The drill apical step = 1.8mmD

Extra Long Drills 16-25mmL

L		Ø	Description	Catalog No.
		1.8/2.5 mmD	Extra-long pilot Drill	DL2525
		2.5/3.2 mmD	Extra-long twist Step Drill	DL2532
		3.2/3.6 mmD	Extra-long twist Step Drill	DL2536

