



UniFit™

Precision you can feel.
Performance you
can trust.

One Unified
Conical Connection
for all Implant Diameters.





UniFit™

Conical Confidence. Unified Platform.

Designed for those who expect more from every implant, every time.

For over 25 years, Adin has been developing implant solutions designed to help clinicians work with greater simplicity, consistency, and predictability. UniFit™ is the evolution of that experience - bringing decades of experience into a smarter, more streamlined implant system.

UniFit™ combines two proven strengths into one intelligent connection - a precise 7° friction-fit conical interface together with a single, unified prosthetic platform that simplifies surgical and restorative workflows.

Experience the evolution of simplicity.

Optimized for Immediate Placement

The implant features a progressive thread design transitioning from coronal trapezoidal, bone-condensing threads to sharp V-shaped cutting threads apically. The design is engineered to promote bone engagement and primary stability for immediate placement and immediate loading protocols

Sized for Versatility. Built for Confidence.

A versatile size range - from 6 mm to 18 mm in length and 3.5 mm to 6 mm in diameter - to cover a broad spectrum of clinical scenarios.



True Friction-Fit Confidence

7° Morse Taper conical connection forms a secure friction-fit between implant and abutment - delivering mechanical stability, minimizing micromovement, and protecting against long-term screw loosening.

Star Index for Surgical Precision

Six-position star index ensures controlled abutment placement.

Built-in Platform Switching

Through its conical connection and anatomically shaped components, UniFit™ incorporates built-in platform switching designed to help preserve marginal bone and support long-term soft-tissue stability.


One Connection for All Sizes

UniFit's smart design simplifies your workflow with a single prosthetic platform across all implant diameters. No juggling components. No guesswork. Just one friction-fit solution for all your cases.






Designed for **Precision.**

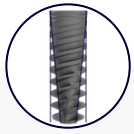

Advanced Implant Design

UniFit's advanced implant design integrates multiple functional zones to support stability across bone types - resulting in enhanced primary stability, bone preservation, and a simplified surgical workflow:

- Trapezoid, bone-condensing coronal threads help reduce crestal pressure and support bone remodeling.
- Sharp, V-shaped apical threads enable efficient bone cutting.
- A tapered core with parallel thread geometry promotes bone engagement while creating healing chambers toward the implant core.
- Self-cutting threads are designed to minimize bone removal and preserve ridge volume.

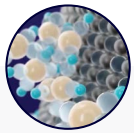


OsseoFix™ Surface Treatment

OsseoFix™ is Adin's advanced calcium-phosphate surface treatment, creating a clean micro-roughened surface designed to support early bone response.


Key Benefits of OsseoFix™:

- Controlled surface roughness for enhanced mechanical interlocking.
- Designed to promote faster woven-to-lamellar bone transition.
- Debris-free finish improves implant-to-bone interface quality.
- Developed to support early-stage stability critical for immediate or early loading.

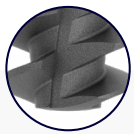

Rounded Apex for Anatomical Safety

The implant's gently rounded tip reduces risk in critical anatomical areas.

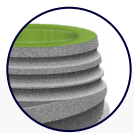

Back-tapered Implant Neck

Designed to reduce crestal bone stress and support long-term marginal bone preservation



Six Position Star Index


A built-in star index enables six stable abutment orientations for controlled placement and improved prosthetic predictability.



Built for Performance.




7° Friction-Fit Conical Connection Precision You Can See at the Interface

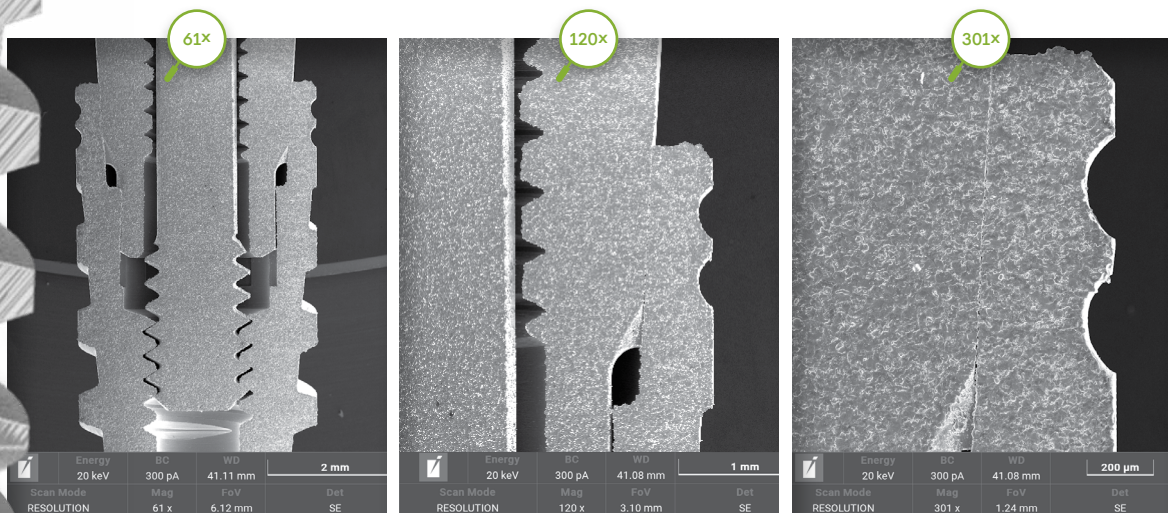

Microscopic view of the implant - abutment interface

- This active connection provides mechanical stability and a tight seal through the taper itself - not just the screw.
- Reduces micromovement and long-term screw loosening.
- Minimizes microgaps that can lead to bacterial infiltration and peri-implant inflammation.

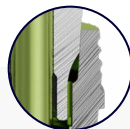

Built-in Platform Switching

Achieved through the conical implant-abutment connection and anatomically shaped restorative components. Research has shown that platform switching helps preserve marginal bone levels and soft tissue stability, supporting better long-term esthetic.



One Connection. Total Restorative Freedom.

UniFit™ is available in diameters from 3.5 mm to 6.0 mm and lengths from 6mm to 18mm, giving you the flexibility to treat a wide range of clinical scenarios – from tight anterior spaces to high-load posterior zones.



	L6mm	L8mm	L10mm	L11.5mm	L13mm	L16mm	L18mm
ø3.5mm		ISPUF0835	ISPUF1035	ISPUF1135	ISPUF1335	ISPUF1635	ISPUF1835
ø3.75mm		ISPUF0838	ISPUF1038	ISPUF1138	ISPUF1338	ISPUF1638	ISPUF1838
ø4.3mm	ISPUF0643	ISPUF0843	ISPUF1043	ISPUF1143	ISPUF1343	ISPUF1643	ISPUF1843
ø5.0mm	ISPUF0649	ISPUF0849	ISPUF1049	ISPUF1149	ISPUF1349	ISPUF1649	ISPUF1849
ø6.0mm	ISPUF0660	ISPUF0860	ISPUF1060	ISPUF1160	ISPUF1360		

Healing Abutments

HealNScan™

	Width/GH	2.0	3.0	4.0	5.0	6.0
	5.0	UF0046	UF0047	UF0048	UF0049	UF0050

Coming
Soon

	Width/GH	2.0	3.0	4.0	5.0	6.0
	Slim (3.5mm)	UF0095	UF0096	UF0097	UF0098	UF0099
	Regular (4.5mm)	UF0090	UF0091	UF0092	UF0093	UF0094
	Wide (5.5mm)	UF0095W	UF0096W	UF0097W	UF0098W	UF0099W

Tools

Implant Drivers					Star Drivers for Prosthetic Attachment					
UF0068	UF0071	UF0074	UF0067	UF0066	RP0064	RP0004	RP0005	RP0065	RP0018	RP0019
										

Screw Pick-Up	Retrieval Tool	Straight TMA				Ball Attachment
RP6400	UF3401	RS6197	RS6196	RS6193	RS6191	RS6090
						

GH Gingival Height RH restorative height

Impression Taking

		Short	Long			Short	Long	Implant Analog	
	Open Tray	UF5610	UF5614		Closed Tray	UF5109	UF5113		UF0007

Straight Titanium Abutments

RH 6.5mm				RH 12.5mm		
Width/GH	1mm	2mm	3mm	1mm	2mm	3mm
Regular (4.4mm)	UF0012	UF0011	UF0010	UF1251	UF1252	UF1253
Wide (5.1mm)	UF0012W	UF0011W	UF0010W	UF1251W	UF1252W	UF1253W
Slim (4.0mm)	-	UF0020 (GH 2.3, RH 6.0)		-		
						

Temporary Abutments (GH 2mm, RH 9.5)

Diameter	Engaging	Non Engaging	Plastic Cylinder
Slim (3.5mm)	UF0036	UF0037	UF0038
Regular (4.5mm)	UF0043	UF0044	UF0045
Wide (5.5mm)	UF0043W	UF0044W	UF0045W
			

Angled Titanium Abutments

Angle	15°			25°		35°	
Width/GH	1mm	2mm	3mm	1.5mm	3mm	1.5mm	3mm
Regular (4.4mm)	UF0013	UF0014	UF0015	UF2515	UF2530	UF3515	UF3530
Wide (5.1mm)	UF0013W	UF0014W	UF0015W	UF2515W	UF2530W	UF3515W	UF3530W
							



One Connection. Total Restorative Freedom.

Digital Restoration

Scan Body	7mm H	10mm H	Digital Analog	
	UF0088	UF0089		UF1000

Ti Bases

RH-4.0mm

		Width/GH	1mm	2mm	3mm	4mm
	Engaged	Regular (4.5mm)	TB07411	TB07421	TB07431	TB07441
		Wide (5.5mm)	TB07411W	TB07421W	TB07431W	TB07441W
	Non Engaged	Regular (4.5mm)	TB07412	TB07422	TB07432	TB07442
		Wide (5.5mm)	TB07412W	TB07422W	TB07432W	TB07442W

RH-6.0mm

		Width/GH	1mm	2mm	3mm	4mm
	Engaged	Regular (4.5mm)	TB07611	TB07621	TB07631	TB07641
		Wide (5.5mm)	TB07611W	TB07621W	TB07631W	TB07641W
	Non Engaged	Regular (4.5mm)	TB07612	TB07622	TB07632	TB07642
		Wide (5.5mm)	TB07612W	TB07622W	TB07632W	TB07642W

RH-8.0mm

		Width/GH	1mm	2mm	3mm	4mm
	Engaged	Regular (4.5mm)	TB07811	TB07821	TB07831	TB07841
		Wide (5.5mm)	TB07811W	TB07821W	TB07831W	TB07841W
	Non Engaged	Regular (4.5mm)	TB07812	TB07822	TB07832	TB07842
		Wide (5.5mm)	TB07812W	TB07822W	TB07832W	TB07842W

Full Arch Solutions

TMA (Trans Mucosal Abutment- Multi Unit)

	Width/GH	1mm	2mm	3mm	4mm	5mm
	Straight	UF0032	UF0033	UF0034	UF3729	UF3730
	Angled 17°	UF3732 (1.5mm)	UF3734 (2.5mm)	UF3737 (3.5mm)	UF3739 (4.5mm)	
	Angled 30°	UF3726	UF3727	UF3735	UF3738	

Ball Attachment

	GH	0.5mm	1mm	2mm	3mm	4mm	5mm
	Ball Attachment	UF0051	UF0052	UF0053	UF0054	UF0055	UF0056

Grip Overdenture Attachment System

	GH	1mm	2mm	3mm	4mm	5mm	6mm
	Grip Abutment	UF0101	UF0102	UF0103	UF0104	UF0105	UF0106

Guided Surgery



Guided Surgery Tools Kit UniFit Implants

GSK005

Surgical Kit



UniFit Surgical Kit Premium

UF6301

UniFit Surgical Kit Advanced


UF6322



Guided Surgery Drills Kit Basic

GSK001

For the
full catalog
Scan here

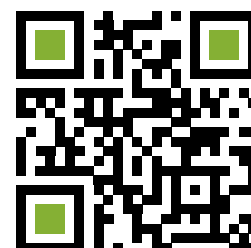



Industrial Park Allon Tavor, POB 1128, Afula 1811101, Israel
T. +972-4-6426732 | www.adin-implants.com